



# Airone FSC-1000/1200PP Forensic Science Safety Cabinet

The Airone FSC-1000/1200PP is a high performance Forensic Science Cabinet designed for the protection of evidence; such as articles of clothing, mobile phones, wallets and bags from particulate contamination and dust.

It provides protection for the operator and the surrounding environment from noxious odours and other contaminants which may emanate from contents held within the cabinet.

The unit comes complete with shower services incorporating a handset with trigger handle, waste outlet and anti-siphon bottle trap. A horizontal hanging rail is also fitted as standard. Optional services are available if required.



## Benefits


- Lockable door
- Simple control panel
- Quiet operation/running mode
- Removable perforated worktop trays
- Shower handset/washdown facility
- Hanging rail
- Visual indication/alarm of airflow
- Hour counter to display running time
- Carbon/HEPA filter technology

Optional internal lighting available

The Airone FSC-1000PP  
1M wide (external width)

The Airone FSC-1200PP  
1.2m wide (external width)

## Technical Specification

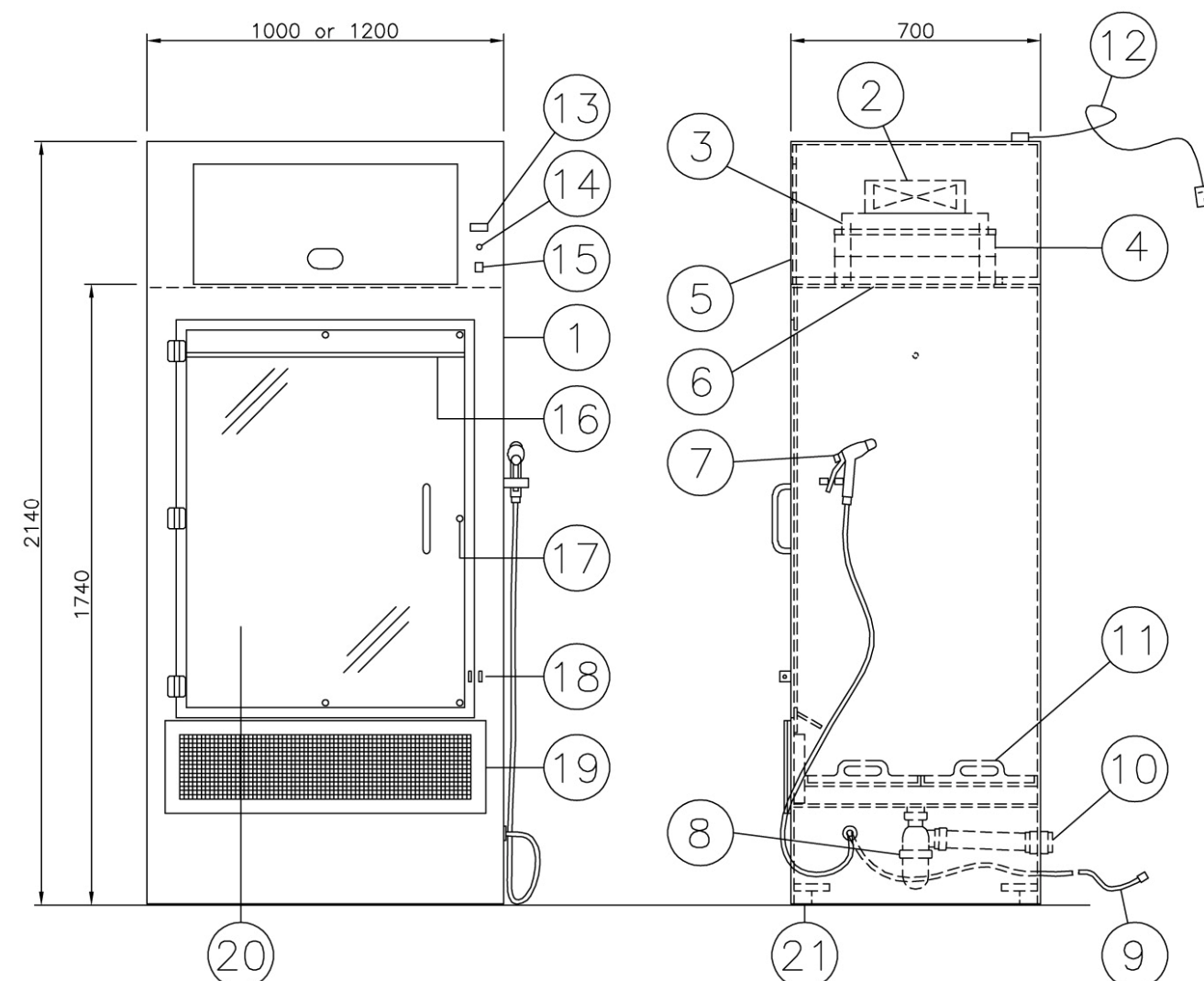

- 1) White polypropylene body
- 2) Exhaust fan
- 3) Filter clamp
- 4) Filters: HEPA & carbon
- 5) Removable access panel
- 6) Filtrete pre-filter
- 7) Shower handset with trigger on handle and flexible shower hose
- 8) Vulcathene glass bottle trap
- 9) Flexible hose with pushfit 1/2" Hozelock connection
- 10) 38mm Vulcathene bulkhead waste fitting

- 11) Removable perforated trays
- 12) Mains lead
- 13) Hours counter
- 14) Low airflow warning LED
- 15) Illuminated green ON/OFF switch
- 16) Hanging rail
- 17) Door catches (1 lockable)
- 18) Security tag location lugs
- 19) Filter
- 20) Clear PVC door
- 21) Adjustable feet

Voltage: 230V 50Hz 1ph

### Airone FSC-1000PP

External dimensions (w x d x h) 1000 x 700 x 2140mm  
Internal dimensions (w x d x h) 975 x 670 x 1380mm

### Airone FSC-1200PP

Ext dims (wdh) 1200 x 700 x 2140mm  
Int dims (wdh) 1175 x 670 x 1380mm