



Recirculates clean air back into room




Airone FSC-1000/1200PP Forensic Science Safety Cabinet

Technical Specification




The Airone FSC-1000/1200PP is a high performance Forensic Science Cabinet designed for the protection of evidence such as articles of clothing, mobile phones, wallets and bags from particulate contamination and dust.

It provides protection for the operator and the surrounding environment from noxious odours and other contaminants which may emanate from contents held within the cabinet.

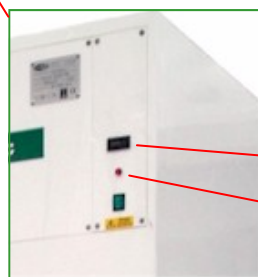


Benefits

- Simple control panel
- Hour counter to display running time
- Visual indication/alarm of airflow
- Quiet operation/running mode
- Two speed fan
- Hanging rail
- Two Shelves
- Removable perforated worktop trays
- White board: for evidence identification
- Shower handset/washdown facility

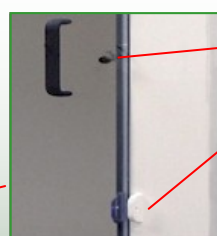

- Lockable door
- Evidence lock
- Pre-filter
- Hose isolation valve (see 9 on opposite page)

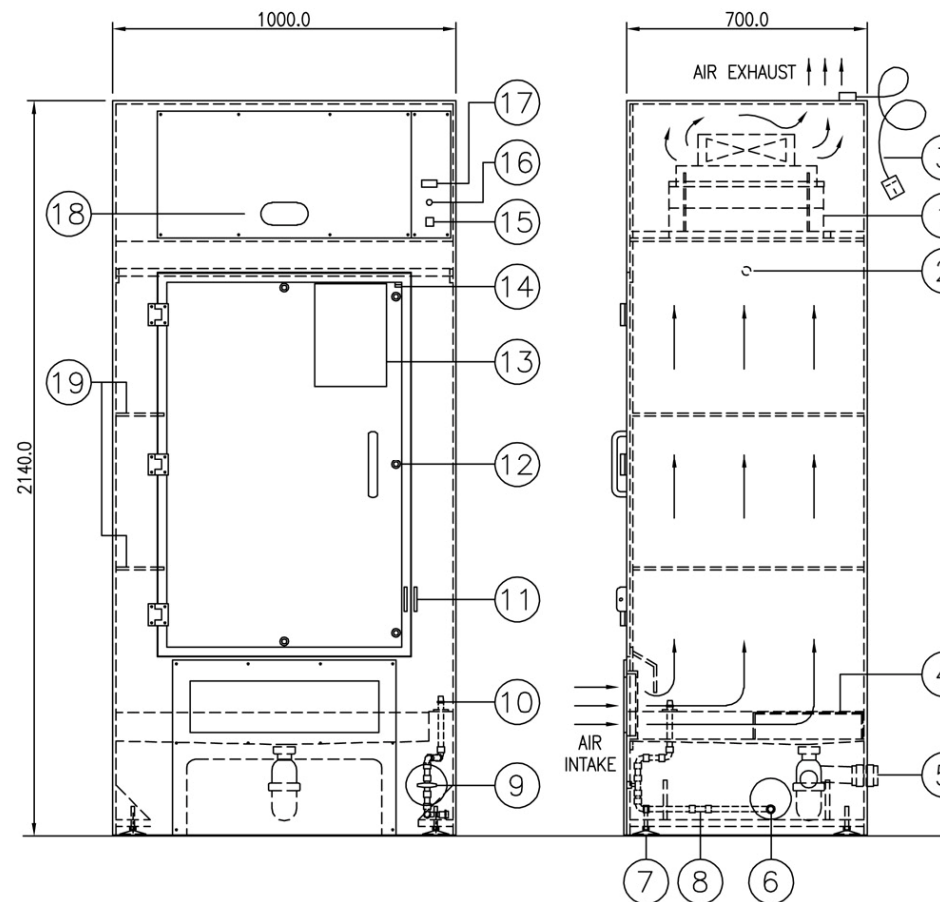

- 1) Carbon and HEPA exhaust filters
- 2) Hanging rail
- 3) Mains Lead
- 4) Removable perforated stand
- 5) Bottle trap & waste outlet (can exit on left, right, or rear)
- 6) Male 3/4" BSP fitting for water inlet
- 7) Adjustable feet
- 8) Water double check valve
- 9) Water isolation valve
- 10) Internal fitting for hose
- 11) Security tag location lugs
- 12) Door catches (1 lockable)
- 13) White board
- 14) Door operated fan speed switch
- 15) Illuminated green ON/OFF switch
- 16) Low airflow warning LED
- 17) Hour counter
- 18) Filter access panel/viewing window
- 19) Perforated shelves

Airone FSC-1000PP

External dimensions (w x d x h) 1000 x 700 x 2140mm
Internal dimensions (w x d x h) 975 x 670 x 1380mm

Airone FSC-1200PP

External dimensions (w x d x h) 1200 x 700 x 2140mm
Internal dimensions (w x d x h) 1175 x 670 x 1380mm

Contact Safelab:

Unit 29 Lynx Crescent
Weston Super Mare
North Somerset
BS24 9BP

Tel: 01934 421 340
E-mail: sales@safelab.co.uk

Fax: 01934 641 569
www.safelab.co.uk

Options

- Internal lighting
- HEPA filter on inlet
- Additional hanging storage

